


# Plan Parameters Report

<b>Hospital Name:</b>	Hospital de Clínicas Caracas
<b>Patient Name:</b>	Duque Lovisone, Laura Mercedes
<b>Patient ID / ID2 / SSN:</b>	12707494 / 16-233 /
<b>Date of Birth:</b>	viernes, 17 de diciembre de 1976
<b>Course ID / Plan ID:</b>	Curso1 / SRS_17_50%
<b>Plan Created:</b>	Jugleys, jueves, 26 de mayo de 2016 11:42 a.m.
<b>Plan Last Modified:</b>	VisionDaemon100, jueves, 26 de mayo de 2016 01:23 p.m.
<b>Plan Approval:</b>	Treatment Approved, jueves, 26 de mayo de 2016 12:56 p.m.
<b>Image ID:</b>	CT_SRS_VMS (image has been processed with Head Frame Detection)
<b>Treatment Unit:</b>	Trilogy
<b>Energy / Dose Rate:</b>	6X-SRS / 1000 MU/min
<b>Fixation:</b>	VMS-Frame
<b>Software Version:</b>	Cone Planning 10.0.42.24719

## Isocenter 1 - Transversal

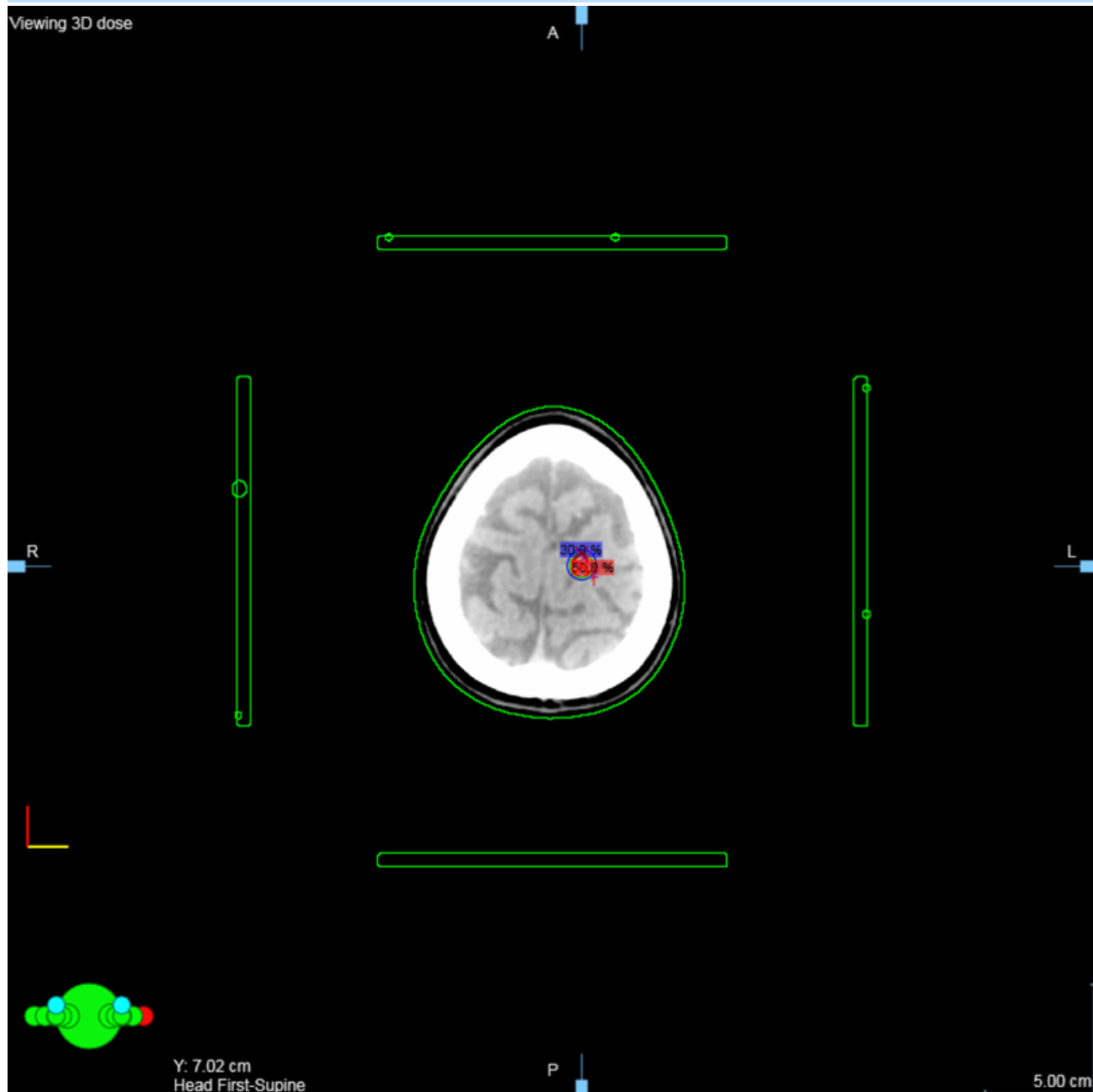

## Calculation

**Calculation Model:** CDC\_10.0.28.2  
**Calculation Options:** Grid: Sparse  
**Arc Angle Resolution:** 10°  
**Dose Grid Resolution:** 1.0mm  
**Calculation Log:** [SRS]  
 Information: Service: ECDC.SRSDose.10.0.28.9833  
 Information: Servant: p2892@e5620\*  
 Information: Client workstation E5620/3728  
 Information: Client version 10.0.42.24719  
 Information: Client time 2016-05-26T11:45:37-04:30  
 Information: Eclipse Cone Dose Calculation (Version 10.0.28)  
 Information: Using DCF protocol version 0.3  
 Information: Beam data directory:  
 \\VARIANDB.ARIA.LOCAL\DCF\$\client\BeamData\CDC\_10.0.28.2  
 Information: Calculated jueves, 26 de mayo de 2016 11:56:26 a.m.

## Prescription

**Dose / Fraction:** 1700.000 cGy  
**Number of Fractions:** 1  
**Total Dose:** 1700.000 cGy  
**Treatment Percentage:** 50.0% of maximum  
**Plan Normalization Method:** No plan normalization  
**Plan Normalization Factor:** 1.000  
**Total Weight:** 15.000  
**Total Rotation:** 1500°  
**3D Dose Maximum (Dmax):** 3400.000 cGy  
**Repeat Factor:** 466.325 cGy

## Parameters

NOTE:

- All angles and field sizes are given using Varian IEC scale.
- Isocenter coordinates are given using the Planning coordinate system.
- Calibration Factors have been corrected to take the absolute dosimetry measurement geometry into account.
- Reference Dose = Repeat Factor × Weight Factor / Average TMR
- Monitor Units = Reference Dose / Calibration Factor

## Treatment Fields

Field ID	Isocenter				Cone [mm]	Couch Rotation	Gantry		Weight Factor	Cal.Fact. [cGy/MU]	Aver. TMR	Ref.Dose [cGy]	Monitor Units	MU/Deg [MU/°]	Aver.d [mm]
	ID	X [cm]	Y [cm]	Z [cm]			Start	Stop							
Campo 2	1	1.37	7.02	-0.72	10	55.0°	20.0°	120.0°	1.000	0.846	0.921	506.383	599 MU	5.99	33.6
Campo 1	1				10	20.0°	20.0°	120.0°	1.000	0.846	0.887	525.679	621 MU	6.21	40.9
Campo 3	1				10	340.0°	240.0°	340.0°	1.000	0.846	0.809	576.522	681 MU	6.81	58.1
Campo 4	1				10	305.0°	240.0°	340.0°	1.000	0.846	0.889	524.698	620 MU	6.20	40.4
Campo 5	1				10	270.0°	240.0°	340.0°	1.000	0.846	0.920	506.908	599 MU	5.99	34.0
Campo 10	2	1.47	6.46	-0.32	5	270.0°	240.0°	340.0°	1.000	0.715	0.879	530.574	742 MU	7.42	38.7
Campo 9	2				5	305.0°	240.0°	340.0°	1.000	0.715	0.842	553.587	774 MU	7.74	45.9
Campo 8	2				5	340.0°	240.0°	340.0°	1.000	0.715	0.765	609.272	852 MU	8.52	63.8
Campo 6	2				5	20.0°	20.0°	120.0°	1.000	0.715	0.853	546.467	764 MU	7.64	44.1
Campo 7	2				5	55.0°	20.0°	120.0°	1.000	0.715	0.885	527.079	737 MU	7.37	37.6
Campo 12	3	1.37	7.63	-0.09	5	55.0°	20.0°	120.0°	1.000	0.715	0.943	494.274	691 MU	6.91	26.3
Campo 11	3				5	20.0°	20.0°	120.0°	1.000	0.715	0.908	513.827	719 MU	7.19	33.1
Campo 13	3				5	340.0°	240.0°	340.0°	1.000	0.715	0.831	561.457	785 MU	7.85	48.3
Campo 14	3				5	305.0°	240.0°	340.0°	1.000	0.715	0.914	510.323	714 MU	7.14	31.7
Campo 15	3				5	270.0°	240.0°	340.0°	1.000	0.715	0.944	493.994	691 MU	6.91	26.3

Cal.Fact. = Calibration factor.  
Aver.d = Average depth.

## Setup Fields

Field ID	Isocenter				Coll Rotation	Couch Rotation	Gantry		Field Size	
	ID	X [cm]	Y [cm]	Z [cm]			Start	Stop	Field X [cm]	Field Y [cm]

Duque Lovisone, Laura Mercedes (12707494)								SRS_17_50% (Treatment Approved)	
G_270	1	1.37	7.02	-0.72	0.0°	0.0°	270.0°	20.0	20.0
G_0	1				0.0°	0.0°	0.0°	20.0	20.0